



RETAJ E



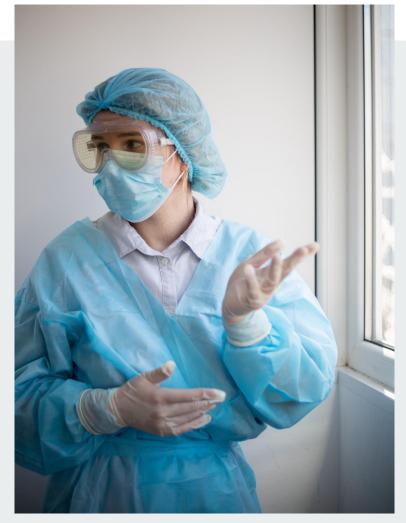

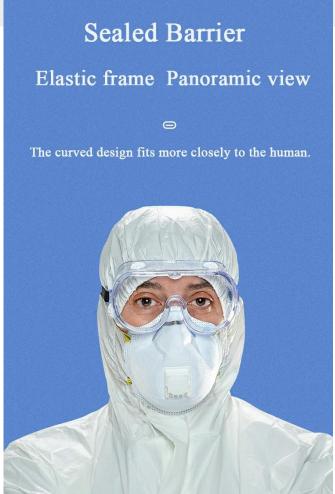

With four holes for ventilation, say goodbye to the foggy surface.

Filter gas and purify air.



# Anti-fog layer lens Not easy to fog

The lens contains an anti-fog layer to prevent fogging of the glasses due to breathing after wearing



# **Products Description**

**Product Type: Safety Glasses Goggles** 

Material: PVC

**Color: Transparent** 

**Size: 18.5\*7.5\*2CM to Variable** 

Weight: 40, 50 grams Certificates: CE / FDA

# Fully enclosed design for wearing with myopia

 $\equiv$ 

Wind and sand dust, blocking blue light and ultraviolet rays.







No burden to wear

Anti-fog material



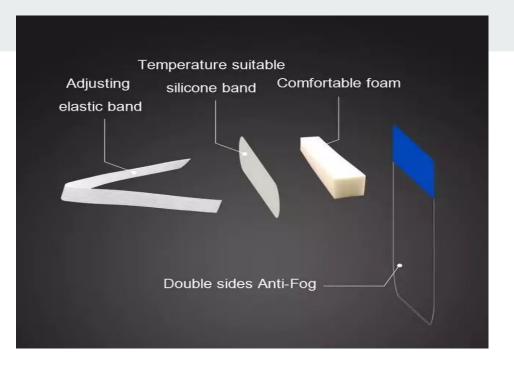

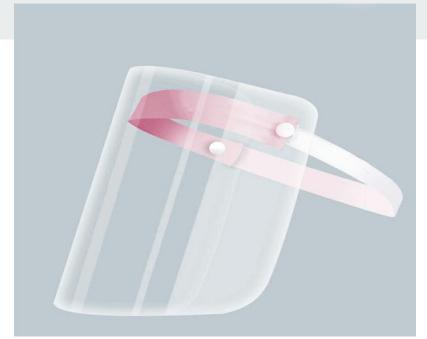

**Products Description** 

**Product Type: Protective Face Shield** 


**Material: Eco Friendly PET** 

Size: 20\*25 cm to Variable

Color: clear

Usage: protect eye,face,mouth Feature: transparent, reusable,

disposable

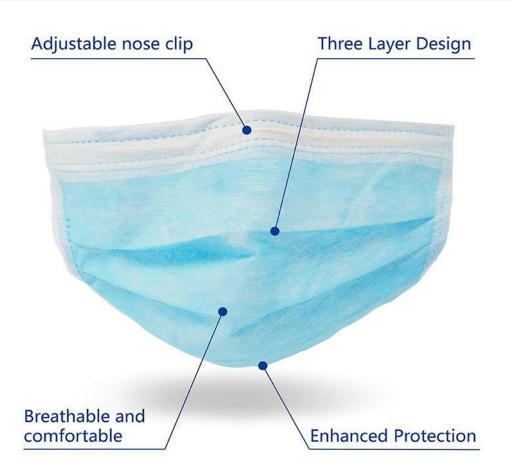

## **Product Description**

**Product Type: 3 Ply Surgical Mask: Sterile and** 

Non- Sterile

Materials: PP nonwoven & filter paper & ES cloth

Size :17.5 x 9.5cm or any you need


Weight: 20g + 20g + 25g or any you need

Color: white, blue, green, pink, etc

Protection: BFE up to 99%, PFE up to 99%

Style : Earloop , ties up , face mask with

shield,face mask spunbond Certificate : CE, ISO ,FDA

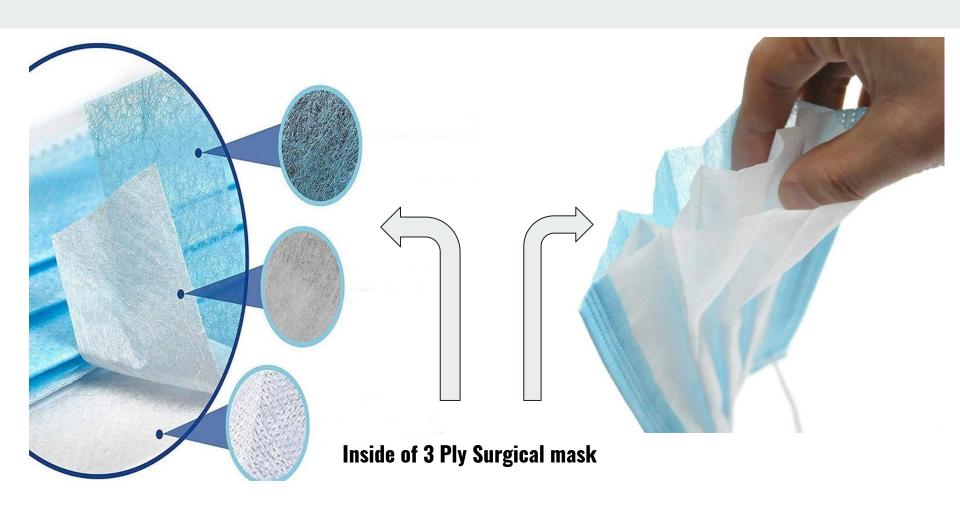

# **N95 MASK**

# **Product Description**

Product: N95 Mask Available Colour: White Function: Dustproof Feature: Comfortable Fit

Type: Disposable

Standard: NIOSH and Non-NIOSH


**Style: Moulded Cup, Cone** 


Material: Non -Woven, Needle Punched Cotton








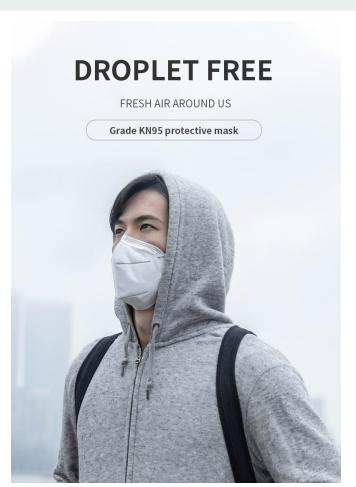


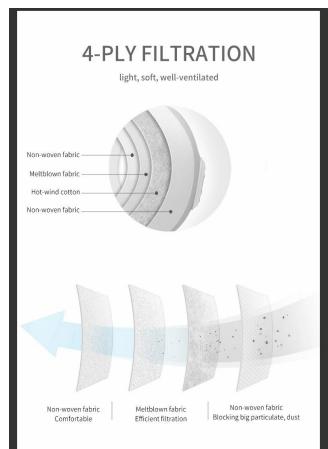


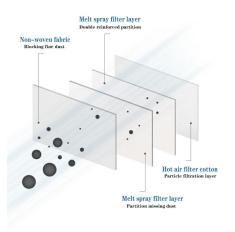





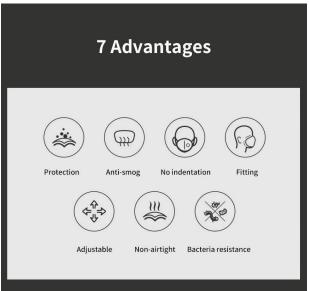

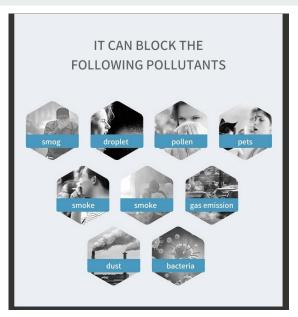

**Product Description** 


Product :KN95 face mask
Material:Polypropylene, Melt blown
Color:White

Feature:Eco-Friendly, Safety, comfortable

Application:Personal protective Style: Elastic Ear Loop or with ties on,Moulded Cup, Cone



#### **WEARING A METHOD**

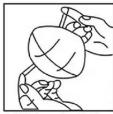

1. Hold the mask so that the stiff bendable strip is on top



2. Pull the bottom of the face mask over your chin.



3. Hold the mask by the ear loops, and place a loop around each ear.



4. Pinch the nose strip to conform to the shape of your nose.



- 5. Check mask tightness before entering work area
- · Cover your face with both hands and exhale quickly
- . If air leaks from the nose, go back to step 4
- . If air leaks from the edges, adjust the earband

Certification/ Class (Standard)	N95 (NIOSH-42C FR84)	FFP2 (EN 149-2001)	KN95 (GB2626-20 06)
Filter performance – (must be ≥ X% efficient)	≥ 95%	≥ 94%	≥ 95%
Test agent	NaCl	NaCl and paraffin oil	NaCl
Flow rate	85 L/min	95 L/min	85 L/min
Total inward leakage (TIL)* – tested on human subjects each performing exercises	N/A	≤ 8% leakage (arithmetic mean)	≤ 8% leakage (arithmetic mean)
Inhalation resistance – max pressure drop	≤ 343 Pa	≤ 70 Pa (at 30 L/min) ≤ 240 Pa (at 95 L/min) ≤ 500 Pa (clogging)	≤ 350 Pa
Flow rate	85 L/min	Varied – see above	85 L/min
Exhalation resistance - max pressure drop	≤ 245 Pa	≤ 300 Pa	≤ 250 Pa
Flow rate	85 L/min	160 L/min	85 L/min
Exhalation valve leakage requirement	Leak rate ≤ 30 mL/min	N/A	Depressurizatio n to 0 Pa ≥ 20 sec
Force applied	-245 Pa	N/A	-1180 Pa
CO <sub>2</sub> clearance	N/A	≤ 1%	≤ 1%

**N95s and KN95s** are both rated to capture 95% of particles, although only KN95 masks are required to pass fit tests. N95 masks have slightly stronger requirements for breathability.















### **Isolation Medical Gown**

Material: PP (polypropylene) Nonwoven, PP+PE material, SMS

material, 25gsm-60gsm

**Cuff: Elastic / Knitted cuffs** 

Collar: with tie or magic stick

Colour: White, Green, Blue, Yellow, Pink, etc.

Sewn seams provide superior strength

Water-resistant, Alcohol-resistant and liquid filter abilities, etc.

Size: S (115x127cm), M (115x137cm), L (120x140cm), XL

(130x150cm)

8. Individual Sterile bag or non-sterile packs are available



## **RETAJ ENTERPRISE LIMITED**

Unit 4 - Unity Works, 329 Petre Street , Sheffield , S4 8LJ

info@retaj.co.uk

+44 7882022903